

# WinRx Dosette/Blister Pack Label Setup Guide

To print dosette labels with med times you must set up a facility record. Go to Browse/Facility and create a new facility with some med times (see example below).

Update Existing Facility

Facility ID  Name  Phone

Address

Notes  Licence #

Days in cycle  Receipt required  Use name order for MARs  Batch minimum interval

**Med Times**

1  2  3  4  5  6  7  8

9  0  A  B  C  D  E  F

G  H  I  J

PRNs on new MAR page  Available only for dispenser type A Separate strip for own bag items  Medadmin leeway  PRN tracking

Bar code on envelope  N=no emr Y=strip emar C=35 day card emar

**Special Times** Breakfast  Lunch  Dinner  Bed  Before breakfast  Before lunch  Before dinner

Emar type


Previous cycle date  Remote printer

Key Commands: F10 = update Esc = cancel

- Go to Browse/System, F2, and check the "Use medtimes" options and hit F10 to save.
- Edit the patient profile and set the facility field.
- When you fill a prescription for that patient, set the med times in the Admin section of the new Rx screen (e.g. "234" would be lunch, dinner, bed).
- Print the dosette labels (in patient rx profile tag rxs with F6, F9 to print). It will look something like the example below:

Rx#	DIN	Drug Qty	Description	Bkf	Lun	Din	Bed
10218	441643	IBUPROFEN 200MG			X		X
J WIEBE take as needed							
10219	291412	TYLENOL 500MG		X		X	X
J WIEBE take as needed							

When the "Use medtimes" option is unchecked it will print like this:

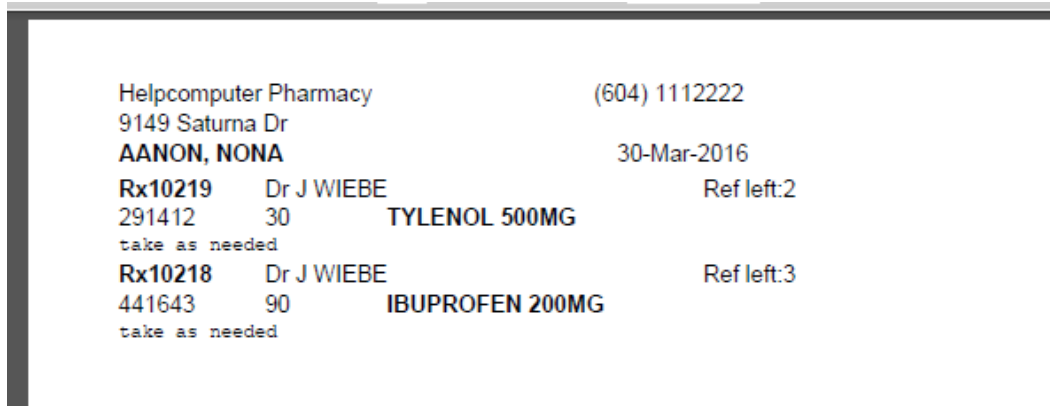

## Compliance Card Setup

To get the compliance card labels that display the quantities in each med time box there are several more steps to set up.

- Go to Browse/Location
- Edit the location and hit F2
- Add a Packaging Files export and set the path to "D"
- You need to have a facility record with all the Special Times fields filled (e.g. 1,2,3,4,5,6,7)
- The patient needs to have a facility set
- The patient needs to have dispenser set to "D"
- Sign in as a supervisor and download the "compliance card" label format in utilities

The prescriptions must have med times set in the "packaging dispense rules" screen for it to work properly (see page 51 of the user guide for more info on that part: <http://arix.ca/files/WinRx-Manual.pdf>).

The printed labels will look something like this:

NONA AANON		Page: 1 of 1		Helpcomputer Pharmacy				
123 ST, BURNABY, BC V4J7A1		31 Mar 2016		9149 Satuma Dr , Burnaby, BC V3J7K1				
604-7779979 BP				604-1112222				
Qty	Drug	Description	Ref Rem	Rx#	Bkf	Lun	Sup	Bed
Instructions								
49	ibuprofen 200mg		3	10218				
take as needed			WIEBE,J	441643	3	2	1	1
35	tylenol 500mg		2	10219				
take as needed			WIEBE,J	291412	2		3	

## Dispill Label Format

Use the same settings as the compliance card setup but select the "dispill" option when installing the label format.